


E-FLEX  **REC**



One file reciprocation for incredible efficiency

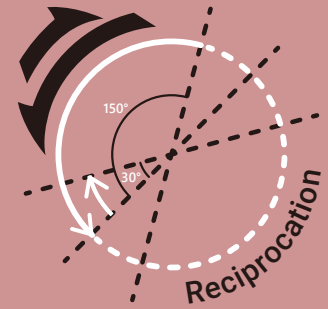

● **ISO Standard Color**

Indicate size of tip

Stopper

Taper Indicator

Precise Length Marking

Variable pitch

Increase gradually from tip to handle to prevent suction and locking



Heat Treated

Control Memory Ni-Ti Wire

Enhance flexibility and resistance to cyclic fatigue for safe treatment

Sharp Cutting Flutes

Enhance cutting efficiency


Double S Cross-Section

Efficient removal of debris

Safety Non-cutting Tip

Reduce ledge formation and apical extrusion, minimize canal transportation

• **Usage Protocol**

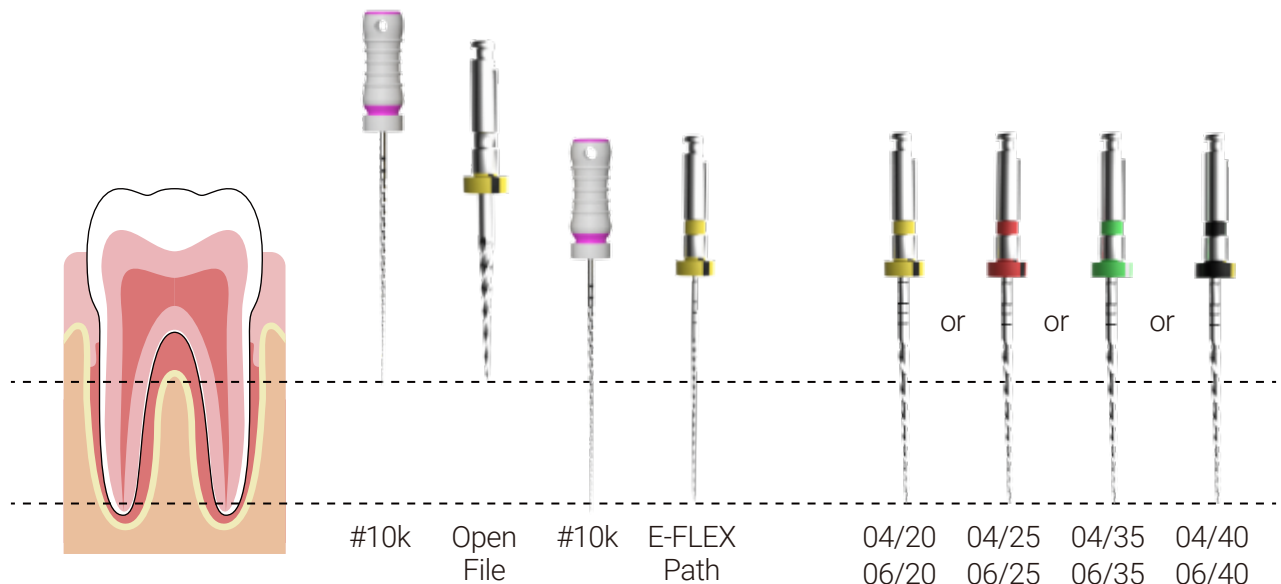



	Color	Tip Size	Taper (%)	Taper Indicator	Cross-Section	Length (mm)	Suggested Angle
E-FLEX REC	●	20	4			21/25/31	RWD: 30° REV: 150°
	●	25	4				
	●	35	4				
	●	40	4				
	●	20	6				
	●	25	6				
	●	35	6				
	●	40	6				