



ENDODONTIC PRODUCTS
CATALOGUE

FANTA
DENTAL

PROFESSIONAL | PRECISELY | RELIABLE

ABOUT FANTA

● Fanta

is an Innovative & High tech manufacturer and designer in the field of Endodontic and Orthodontic material

For different clinical requirements, our product line is so complete . We develop and design our own research including shedding, dredging, preparation, re-treat-ment and files especial for deciduous teeth, for bypass a ledge and for cleaning.

We issued the first preparation kits containing path file in 2014, issued preparation kits containing two sets of path file in 2017, and preparation kits containing three sets of path files in 2019. What's more, we pay more attention for safe dredging and high-efficiency preparation.

Exclusive Agent in Syria **Alkaram Dental Supplies**

 Aleppo-Syria  0212666929  0940666929 Aleppo  0962233550 Lattakia

  @AlkaramDentalSupplies  <http://alkaramdent.com>  0962233622 Damas

PRODUCTS

FANTA

ROTARY & RECIPROCATING FILES

F-ONE	01
F-ONE	02
S4 AF BLUE	03
RETREATMET	04
AF™ BLUE R3	05
C PATH FILES	06
AF™ ROTARY	07

ACCESS & HAND FILES

K FILE	09
H FILE	09
C FILE	09

OBTURATION & MATERIALS

PLUGGER	11
---------------	----

F-ONE

Available sizes:

F one 35/04 , F one 30/06,
F one 30/04 , F one 25/06,
F one 25/04 , F one 20/06, F one 20/04



1. New generation of one file systems
2. Better efficient cutting, time and effort efficiency and easy to use
3. One file is enough to complete the shaping and apply least efforts on the root due to cross sectional of the file
4. AF™-R wire; Improved resistant; 600% higher resistance to cyclic fatigue than normal Ni-Ti wire
5. Inactive Tip: Advanced tip process to avoid forming steps; Suitable for curved canals
6. Magic Flat Design:
Provides space inside the canal to store irrigation fluid and collect debris to extract them to orifice Reduce cutting edges along with the file on one side to reduce cyclic fatigue on the file
7. 11mm file shank that makes F one suitable for molar treatment with limited mouth opening
8. All Fanta files are heat treated and rewarded ISO and FDA



F-ONE

Available sizes:

F one 35/04 , F one 30/06,
F one 30/04 , F one 25/06,
F one 25/04 , F one 20/06, F one 20/04




1. New generation of one file systems
2. Better efficient cutting, time and effort efficiency and easy to use
3. One file is enough to complete the shaping and apply least efforts on the root due to cross sectional of the file
4. AF™-R wire; Improved resistant; 600% higher resistance to cyclic fatigue than normal Ni-Ti wire
5. Inactive Tip: Advanced tip process to avoid forming steps; Suitable for curved canals
6. Magic Flat Design:
Provides space inside the canal to store irrigation fluid and collect debris to extract them to orifice Reduce cutting edges along with the file on one side to reduce cyclic fatigue on the file
7. 11mm file shank that makes F one suitable for molar treatment with limited mouth opening
8. All Fanta files are heat treated and rewarded ISO and FDA

S4

AF BLUE

ROTARY FILES

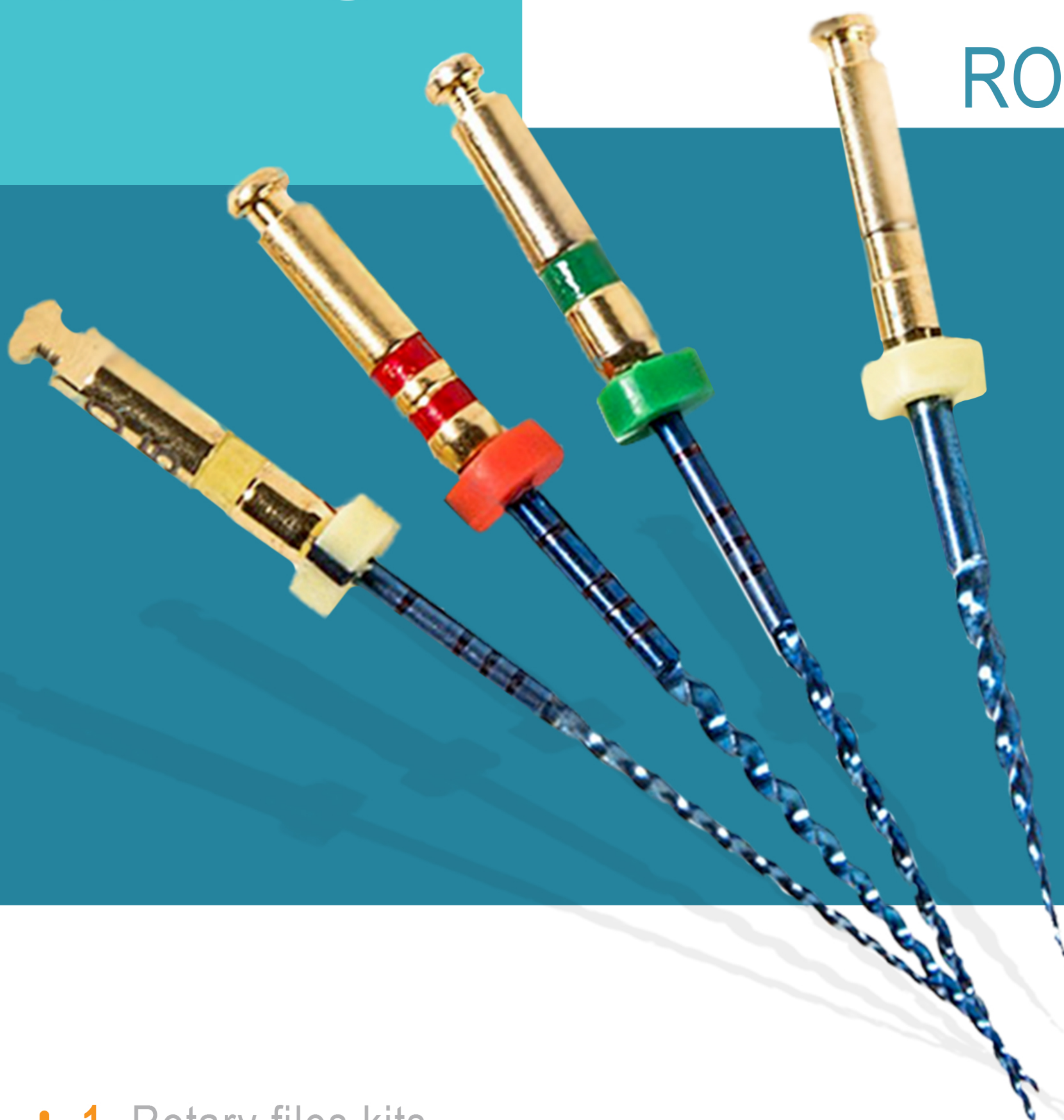


1. Rotary files kits
2. Heat treated files with oxide titanium surface layer which gives the file more resistance and flexibility to prepare narrow canals with multiple file systems, in addition to its double spirals for more cutting efficiency
3. Provides suitable preparation of the canal and it contains the following files:

SX orifice opener	#18 5%
#25 6%	#35 4%

RETREATMET


High efficiency and safety

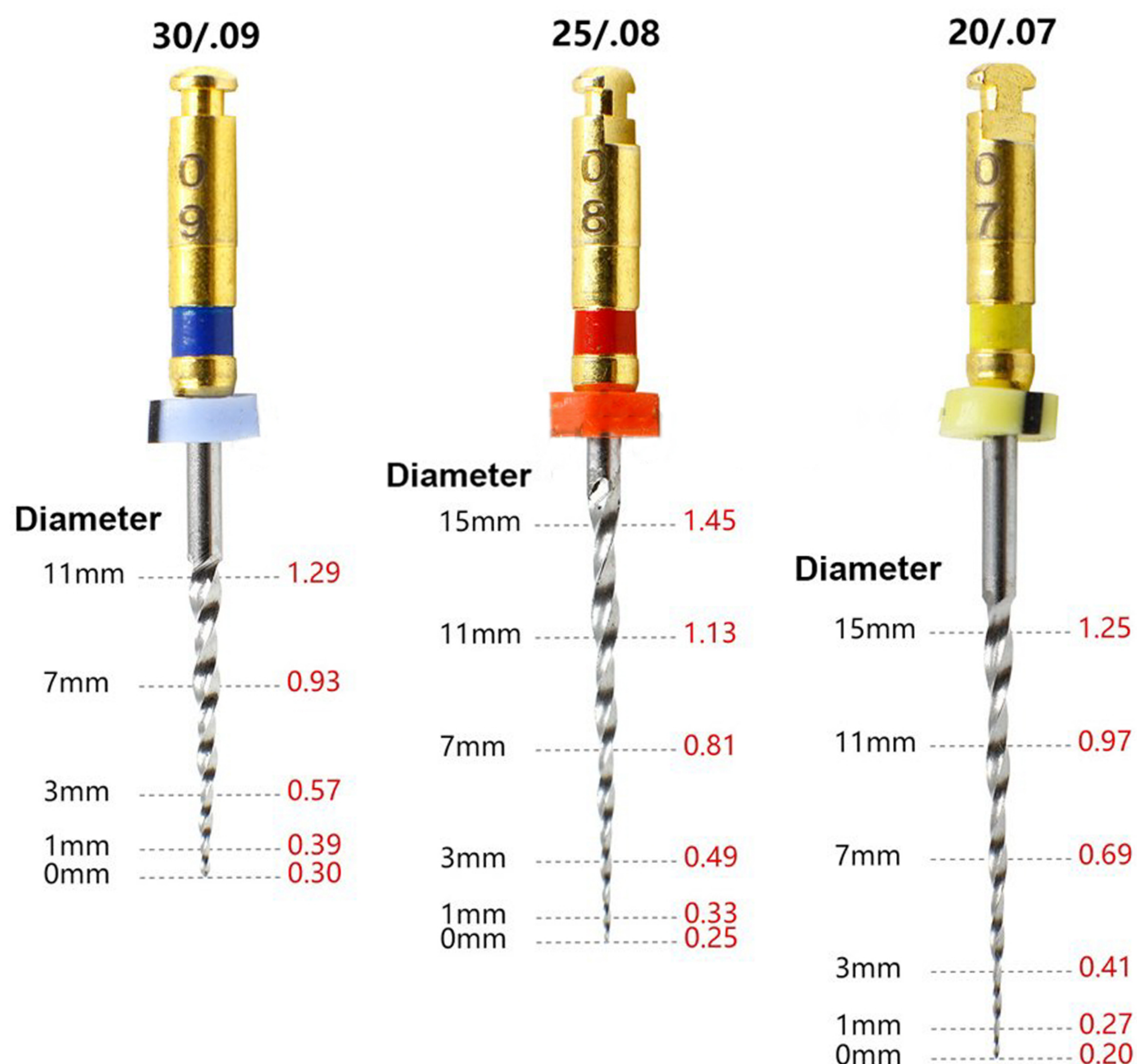

1. New level of root canal retreatment that provide high efficiency and safety in retreatment cases
2. Active Tip system
3. Provides spaces along the file that accumulate the debris & push it out toward the orifice
4. Use file (#30/0.09) 16mm, (#25/0.08) 18mm then file (20/0.07) 22mm
5. Speed: 300-350 rpm torque 2N
6. Advanced control memory alloy tech, AFTM-H wire; Improved resistant; Much better cyclic fatigue resistance

ROTARY FILES

AFTM BLUE R3

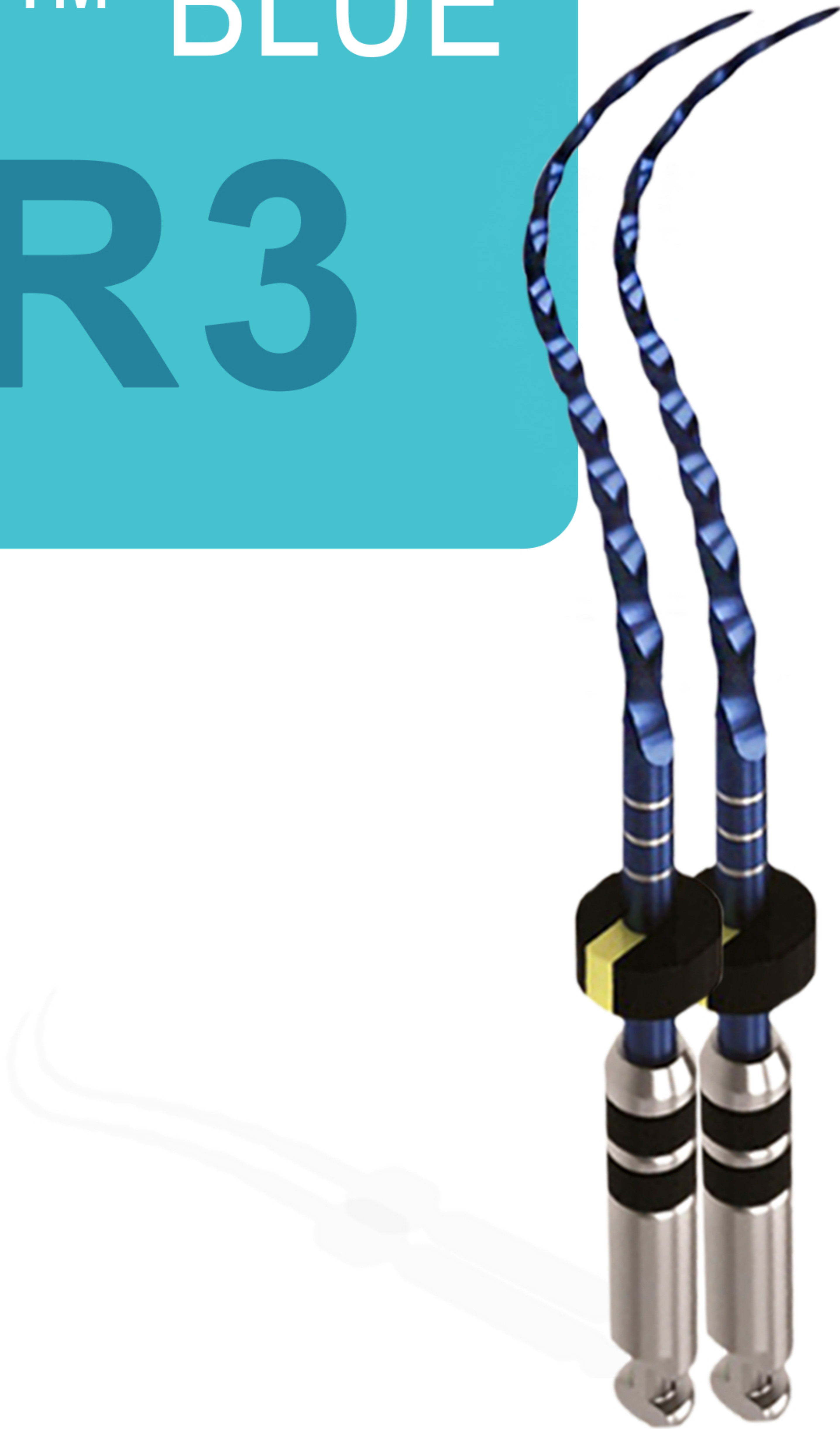

- AFTM-H Wire Tech
- Counterclockwise to 150 degrees, clockwise by 30 degrees
- 360 degrees counterclockwise when used are excellent debris discharge capability
- Increase the file resistance to cyclical fatigue
- AF BLUE Resiproc filed R3 have the same square cross-section and the groove sections upwards

C-PATH ROTARY

ROTARY FILES



Available in 3 colors

- Provides safety negotiation in narrow and curved canals
- AFTM-L wire; Much better Cyclic fatigue resistance, inactive tip
- More flexible than K-FILE
- Speed : 300 RPM torque : 2 N
- Available in 3 sizes:


#13 2%

#16 2%

#19 2%

AFTM RPTARY

Excellent Adaptation to Severely Curved Canals

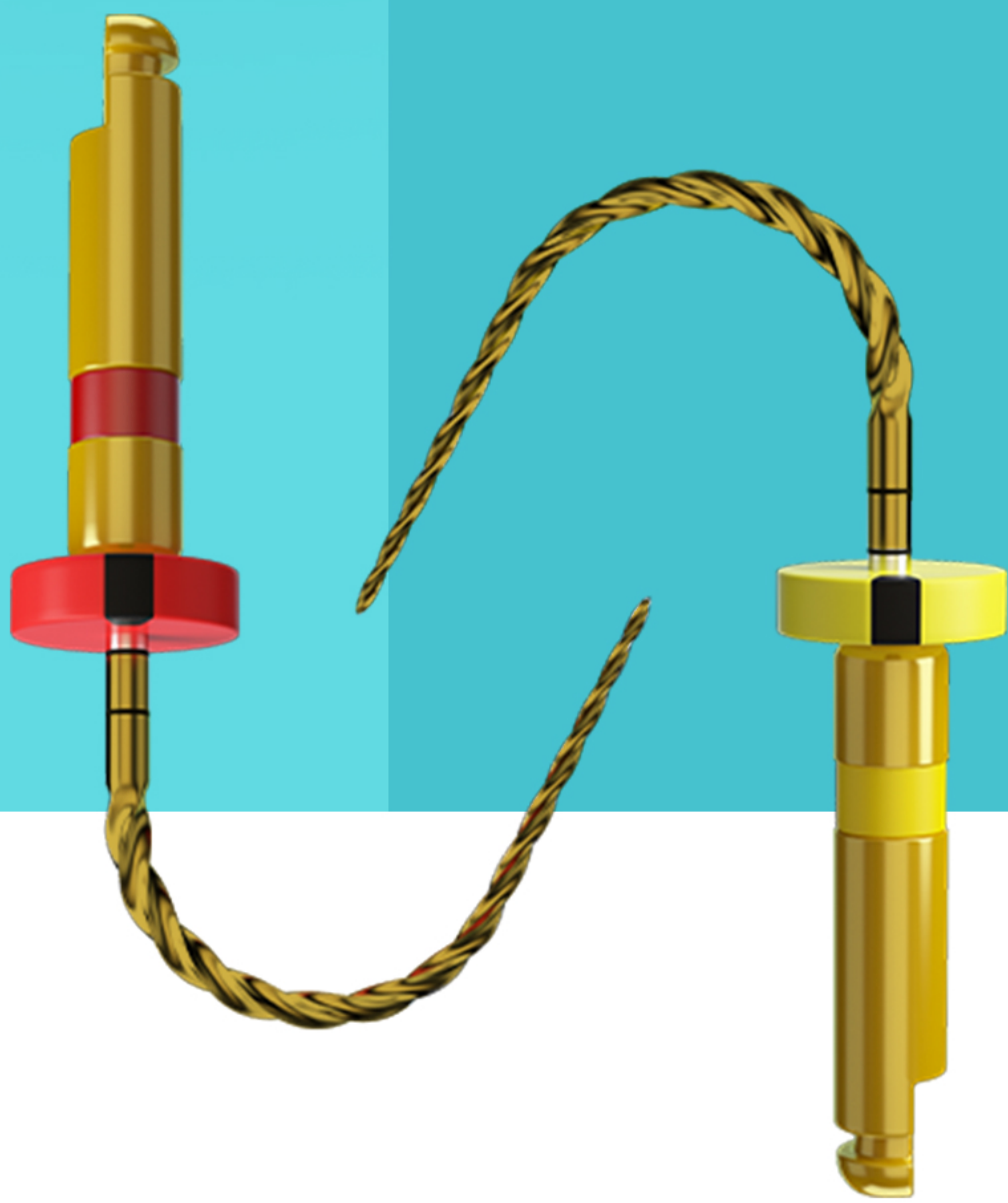

- Progressive Endodontic file system
- AF-H Wire technology
- Extremely flexible files
- Files of choice for very narrow and calcified canals
- Files of choice for severely curved canals
- Minimal invasive files
- Improved resistance to cyclic fatigue





HAND FILES

K FILE | H FILE | C FILE

Hand file is the most widely used ISO standard preparation instrument in clinical practice & has a strong cutting ability. And our files add one more element to take care of the convenience based on the original hand file, rubber handle, we choose the special material for slip prevention better when dentists hold the handle.

K FILE

FOR K FILE

It is a kind of stainless steel wire twisted into a spiral, & its working end is composed of a spiral cutting ridge. The larger the size of K file, the smaller the cutting angle, & the denser the spiral grain, the higher the cutting efficiency.



H FILE

FOR H FILE

It is so similar to the tail of a mouse and named a rat tail, it has a strong cutting ability, so it is generally used for lifting up & down. Stainless steel instruments have advantages over nickel-titanium instruments in relatively straight and irregular root canals.



C FILE

FOR C FILE


It is especially suitable for extremely curved calcified root canals. It increases the possibility of negotiating calcified curved root canals. The success rate of canal treatment was further improved. Due to the unique manufacturing process & the vacuum condition of smelting, the tip has high strength & good flexibility at the back, which greatly reduces the occurrence of root canal transportation or lateral penetration.



OBTURATION & MATERIALS



PLUGGER

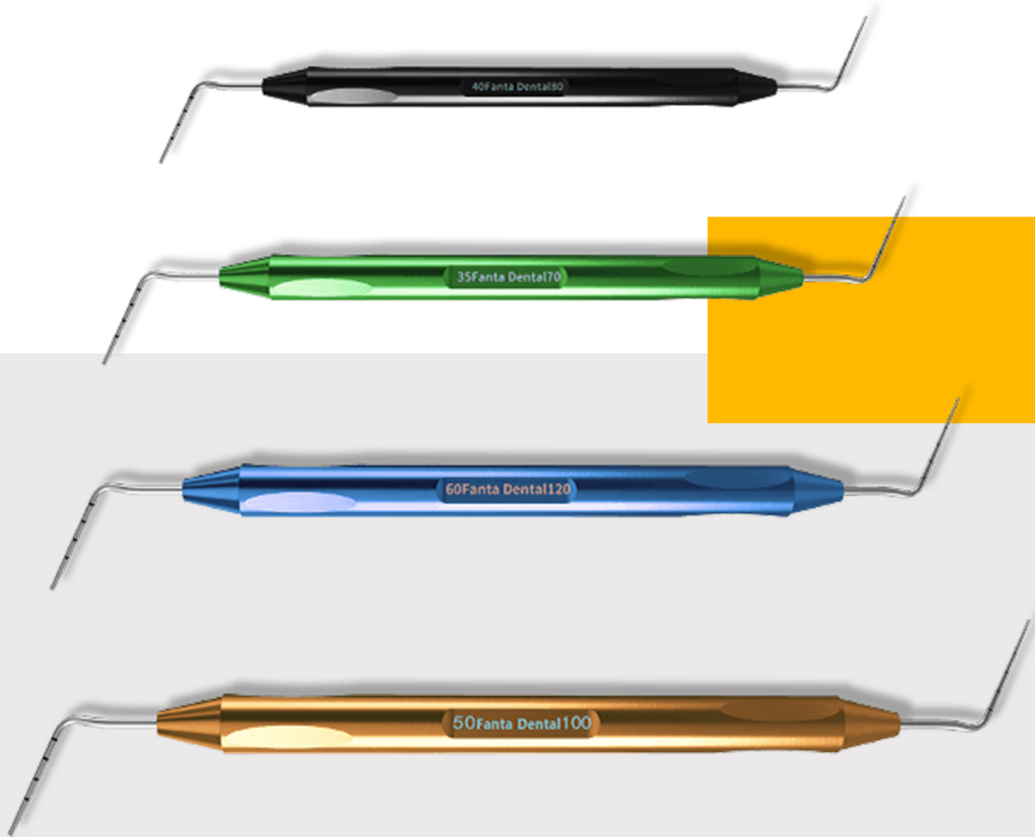

- Stainless steel plugger with length markings
- Flexible Ni-Ti hand plugger for narrow and curved canals
- ISO color marking
- Specialty for Both NITI and Stainless Steel Tips
- Double tip for two different sizes



FANTA DENTAL

Exclusive Agent in Syria **Alkaram Dental Supplies**

🏠 Aleppo-Syria 📞 0212666929 📱 0940666929 Aleppo 📱 0962233550 Lattakia

📱 @AlKaramDentalSupplies 🌐 alkaramdent.com 📱 0962233622 Damas